




doğsan

surgical sutures

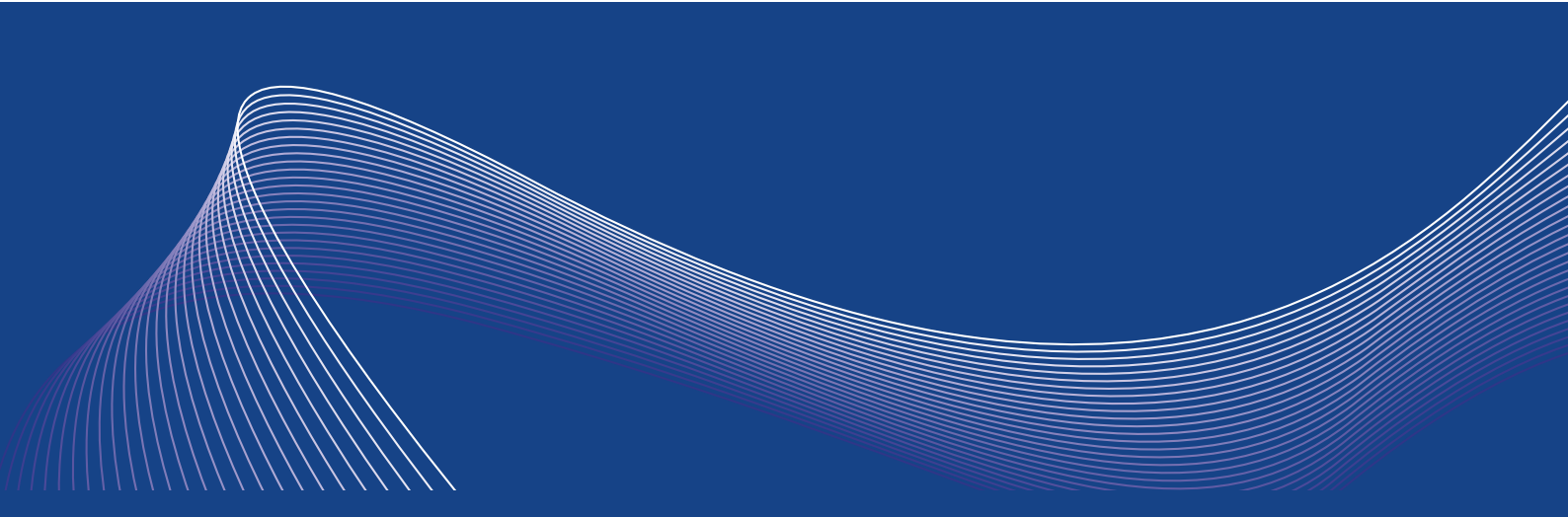

doğsan



Continuous Improvement

The main goal of Doğsan is increasing customer satisfaction through continuous improvements.

Teamwork

Doğsan employs highly experienced and skilled personnel and believes in working as a team to accomplish its goals.

your reliable partner



Experience

Established in 1970, Doğsan has a long tradition in suture manufacturing and is the first sterile surgical suture manufacturer in her region.

Doğsan is a well respected brand name in her region for quality and is now offering its quality products to global markets.

Sutures to meet every surgeons needs

Doğsan has a comprehensive range of absorbable and non-absorbable sutures. Manufacturing sutures ranging from USP 11/0 to USP 8, Doğsan offers an extensive variety of products for all branches of surgery including, cardiovascular surgery, microsurgery and ophthalmology.

Customer Focus

Doğsan is committed to understand and meet the needs of her customers: Surgeons, Nurses and Distributors.

As the requirements of modern surgery change, Doğsan uses its experience and customer focus to customize its technology and product range to meet the most demanding surgeons and nurses.

Commitment to Global Quality Standards

We are committed to global quality principles. Our main goal is to increase customer satisfaction and continuously improve our quality by working as a team.

EP (European Pharmacopoeia)

USP (United States Pharmacopoeia)

ISO 13485:2016

2017/745 Medical Devices Regulation,



your reliable partner

Product Index

absorbable



non - absorbable




microsurgery and
ophthalmology

How To Use This Product Guide


Needle description - Needle type and the description of the needle		Needle profile	Needle diameter	Colour of suture strand	Length of the suture (cm)	Size Grid - Indicates the USP size (below) and metric size (top) for each product listed.		
needle information			needle diameter	colour	suture length	0.7 6/0	1 5/0	1.5 4/0
13mm								
3/8 round bodied		330 μ	undyed	45 cm	J131			
		330 μ	violet	45 cm	L131			
		330 μ	violet	60 cm		L231		
		430 μ	violet	60 cm			L331	
3/8 PL		380 μ	violet	45 cm	L139			
		430 μ	violet	45 cm		L239		
needle length						Product Code		

Explanation of Labelling

Product Label

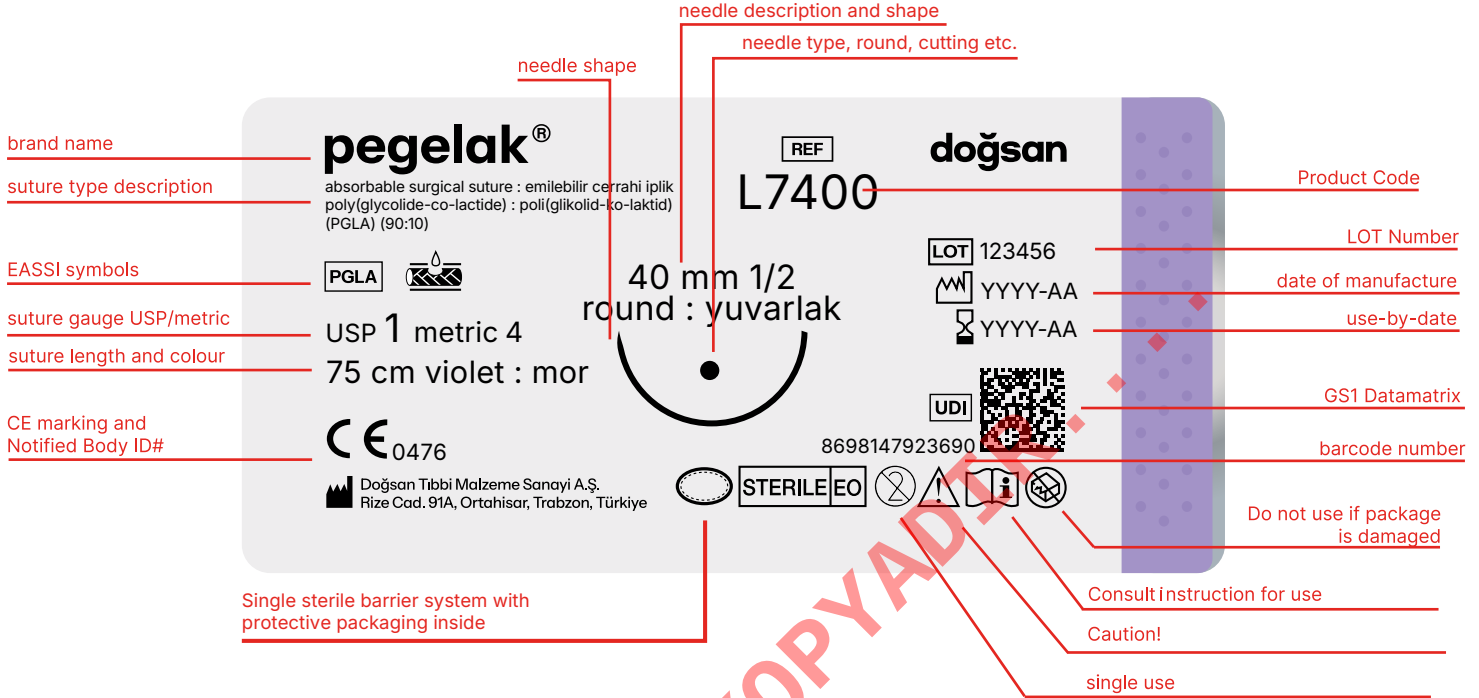

Dozen Box Label

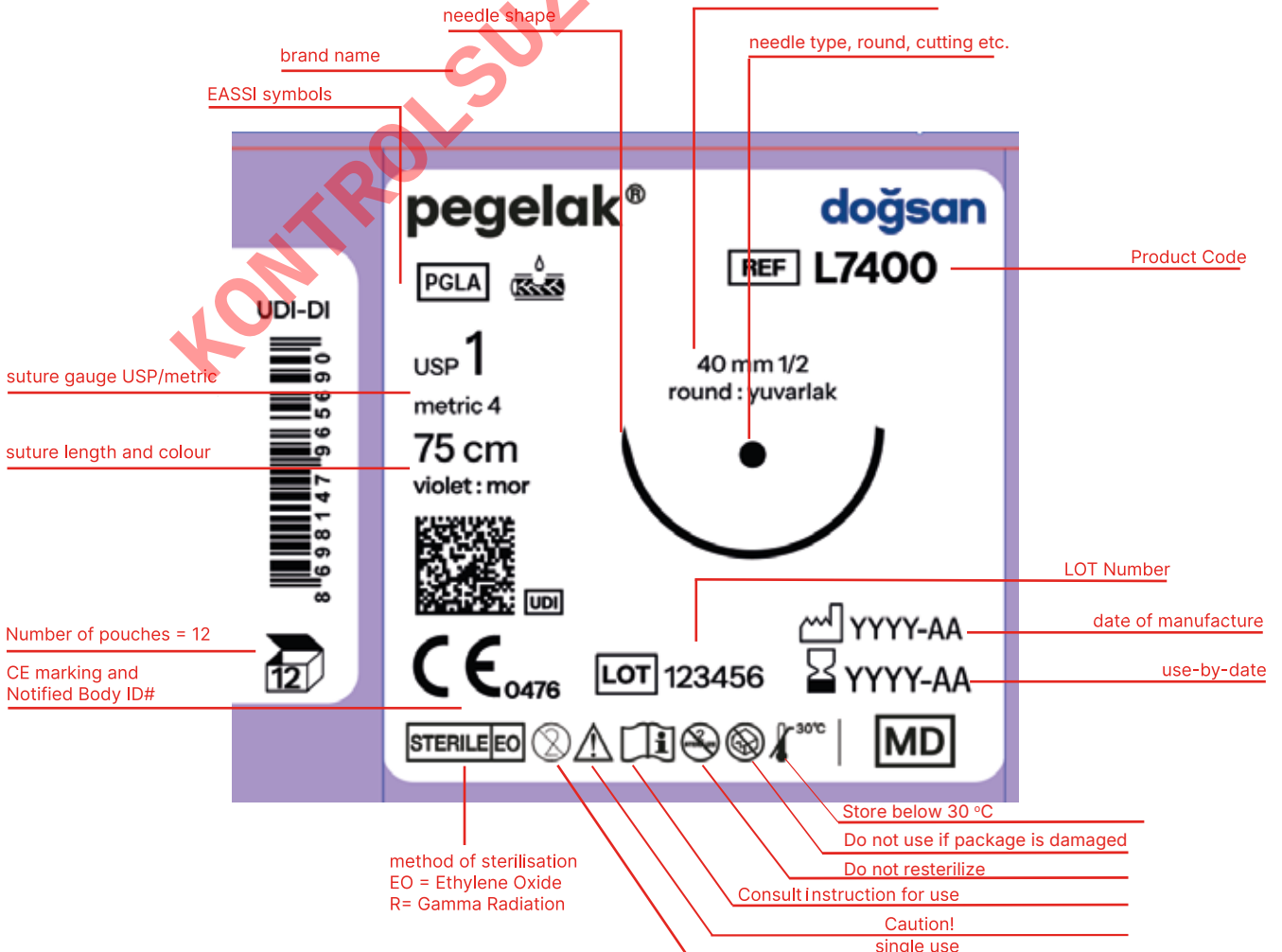


EP/USP Tensile Strength Table (minimum values)

		NON-ABSORBABLE					SYNTHETIC ABSORBABLE			
		DIAMETER	TENSILE STRENGTH		STEEL WIRE		DIAMETER		TENSILE STRENGTH	
		(mm)	(kgf)		(kgf)		MONOFILAMENT	MULTIFILAMENT	(kgf)	
USP	metric	min-max	min	avarage	min	avarage	min-max (mm)	min-max (mm)	min	avarage
11/0	0.1	0.010 - 0.019	-	0.003			-	-	-	
10/0	0.2	0.020 - 0.029	-	0.010			-	-	-	
9/0	0.3	0.030 - 0.039	0.006	0.035			-	-	0.023	0.045
8/0	0.4	0.040 - 0.049	0.015	0.060			-	0.040 - 0.049	0.035	0.070
7/0	0.5	0.050 - 0.069	0.030	0.102			0.050 - 0.094	0.050 - 0.069	0.070	0.140
6/0	0.7	0.070 - 0.099	0.060	0.150			0.095 - 0.149	0.070 - 0.099	0.130	0.250
5/0	1	0.100 - 0.149	0.102	0.310			0.150 - 0.199	0.100 - 0.149	0.340	0.680
4/0	1.5	0.150 - 0.199	0.150	0.510			0.200 - 0.249	0.150 - 0.199	0.470	0.950
3/0	2	0.200 - 0.249	0.310	0.920			0.250 - 0.339	0.200 - 0.249	0.890	1.750
2/0	3	0.300 - 0.349	0.920	1.530			0.340 - 0.399	0.300 - 0.349	1.340	2.680
0	3.5	0.350 - 0.399	1.330	2.240			0.400 - 0.499	0.350 - 0.399	1.850	3.900
1	4	0.400 - 0.499	1.530	2.750			0.500 - 0.570	0.400 - 0.499	2.540	5.080
2	5	0.500 - 0.599	2.240	3.570	4.420	5.900	0.571 - 0.610	0.500 - 0.599	3.180	6.350
3 & 4	6	0.600 - 0.699	2.750	5.100	6.830	9.110	-	-	-	-
5	7	0.700 - 0.799	3.570	6.320	8.550	11.400	-	-	-	-
6	8	0.800 - 0.899	5.100	7.450	10.200	13.600	-	-	-	-
7	9	0.900 - 0.999	-	-	11.930	15.980	-	-	-	-
8	10	1.000 - 1.099	-	-	13.650	18.200	-	-	-	-




Needles




Structure of Atraumatic Needles

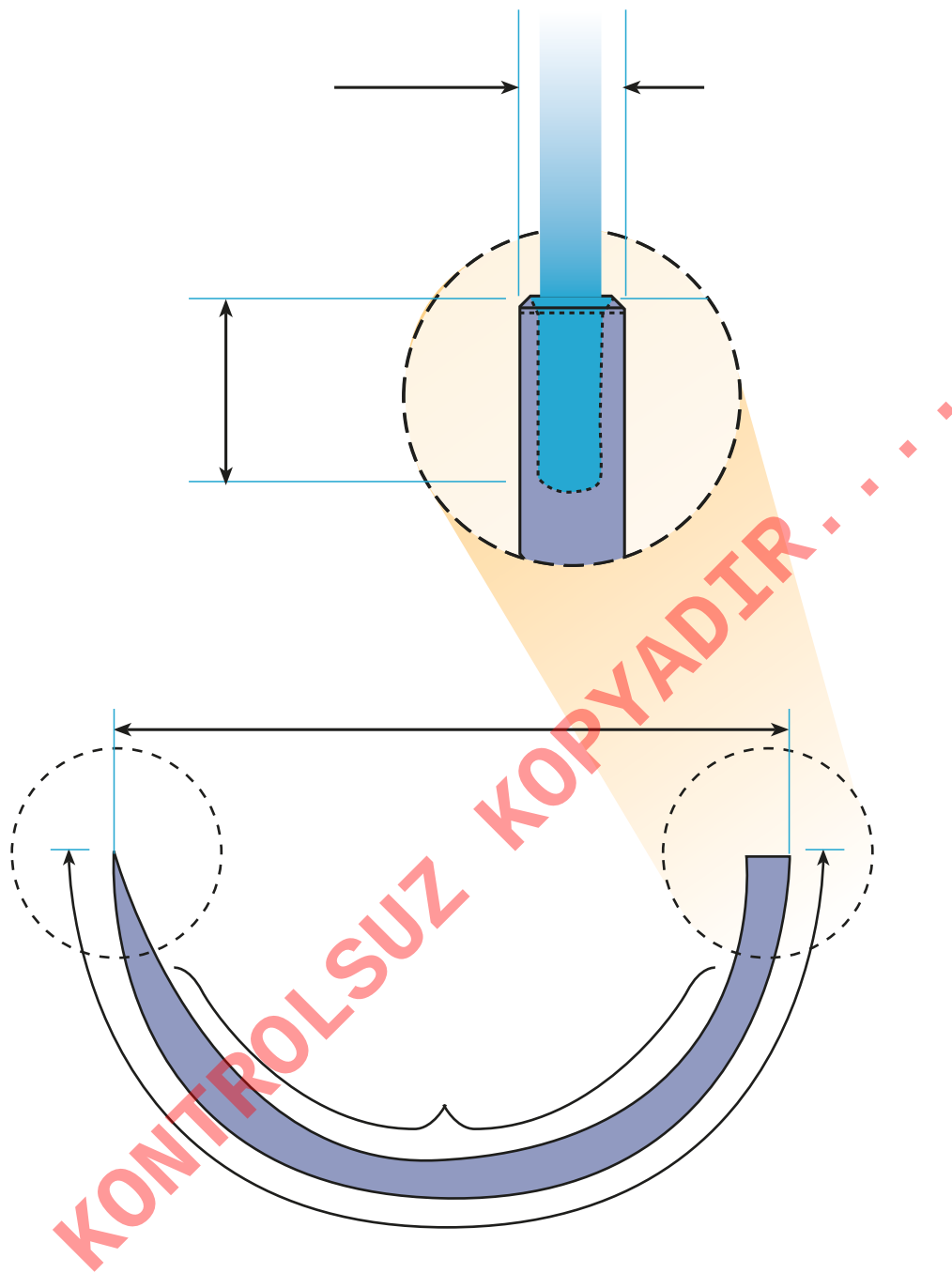

Needle Types According to Curvature


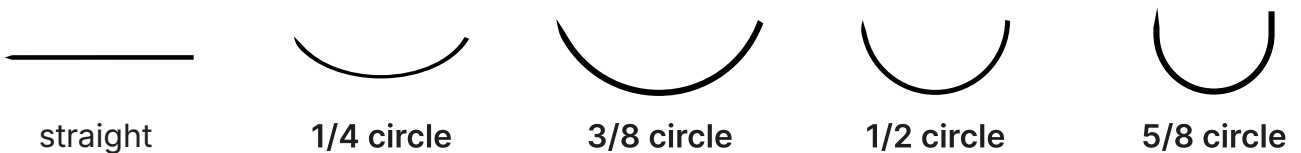


Atraumatic Surgical Needles

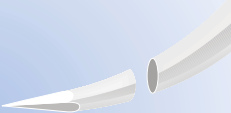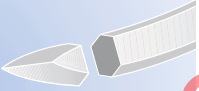
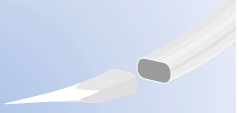
name	tip	body	symbol	description
round bodied				- The most popular and classic needle type. - Use in all soft tissue approximation including cardiovascular surgery.
cutting	inside 			- Edges are specially sharpened for an easy and smooth penetration. - Cutting edge is inside the needle. - Used in skin closures and in all applications requiring cutting needles.
	outside 			- Edges are specially sharpened for an easy and smooth penetration. - Cutting edge is outside the needle. - Widely used in skin closures and especially in plastic surgery.
KD 				- Specially designed for cardiovascular surgery for calcified vessels and tough fibrous tissues. - The body is pressed to square shape to increase resistance to bending. - Tip is specially sharpened for superior penetration.
black 				- Used where better visibility and prevention of reflection of the needle is required (eg when blood concentration is high). - Darkening of the metal is made by surface oxidation technique. - Widely used in cardiovascular surgery.
spatuled				- Both side edges and upper surface are cutting while the body is flattened. - Widely used in ophthalmic surgery.
blunt point				- The tip of the needle sharpness of round bodied is reduced. - Uses include repair of liver and sensitive friable tissues.
diamond (taper cut)				- For better penetration and less tissue trauma, tip of the needle is sharpened. - Widely used in cardiovascular/thoracic, orthopedic, obstetrics/gynecologic surgery and repairs of tendons and fibrous
PL 				- Special geometry for easy penetration and high resistance to bending with a slim body. - Widely used in plastic and cosmetic surgery and skin closures.

Round Bodied (taper) Needles

3/8

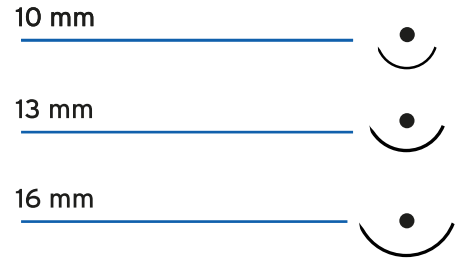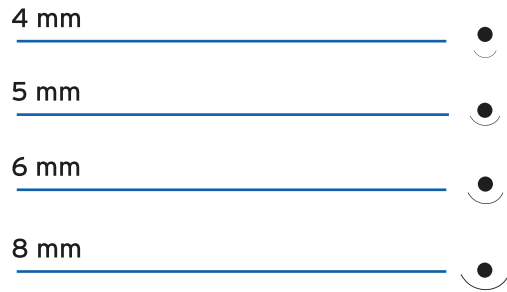

1/2

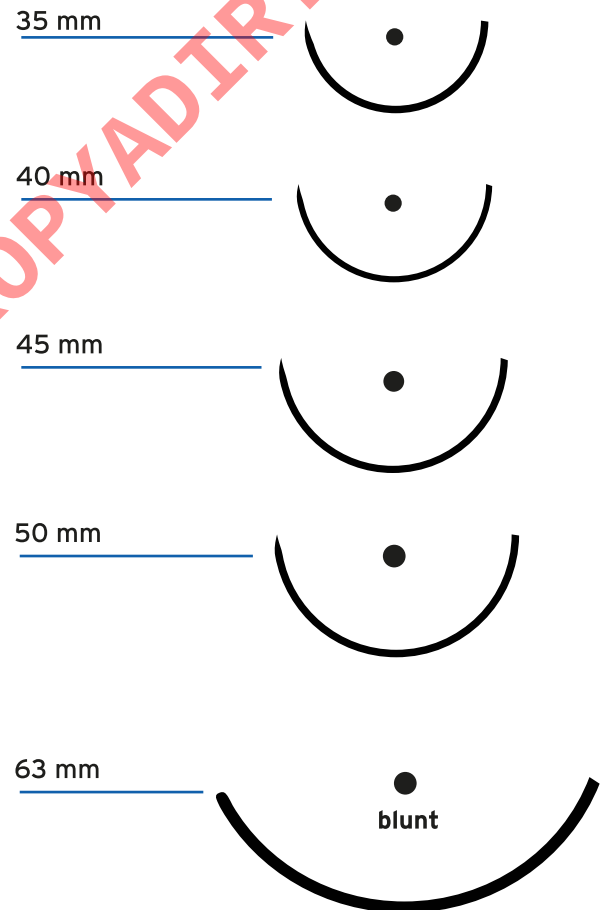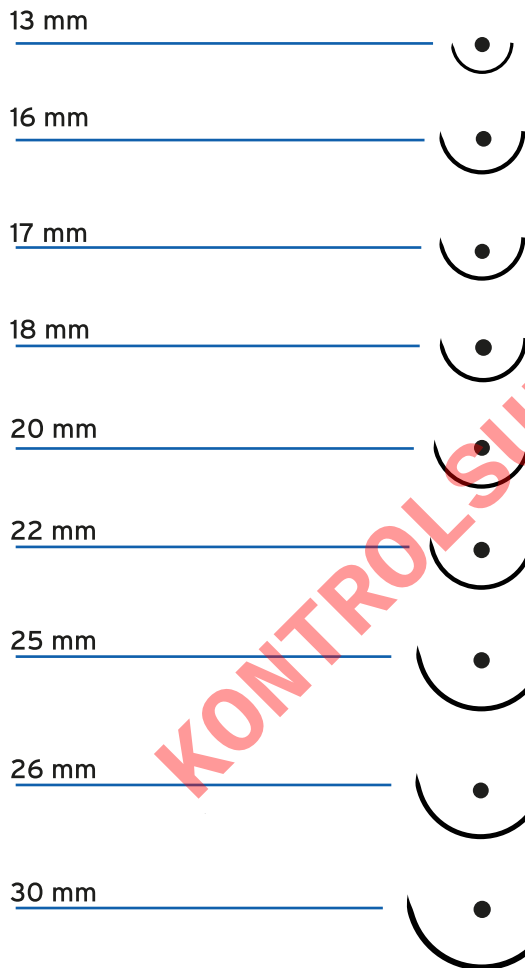

5/8



Special Needles

Diamond
(taper cut)

KD

3/8

18 mm



1/2

17 mm



25 mm



40 mm



45 mm



55 mm



3/8

8 mm



10 mm



10 mm



black

13 mm



1/2

16 mm



20 mm



25 mm



Special Needles

Black

3/8

8 mm 
black

10 mm 
black

10 mm 
black

13 mm 
black

1/2

16 mm 
black

20 mm 
black

25 mm 
black

Spatulated Needles

3/8

6 mm 

1/4
8 mm 

PL Needles

13 mm 

16 mm 

19 mm 

25 mm 

Cutting Needles

1/2


22 mm

30 mm

37 mm

40 mm

50 mm

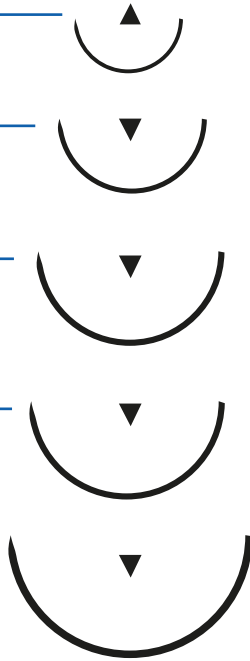


straight

13 mm

18 mm

50 mm

60 mm

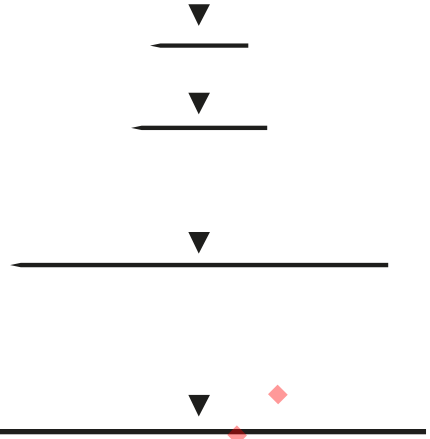

3/8

8 mm

10 mm

13 mm


16 mm

18 mm

20 mm

25 mm

26 mm

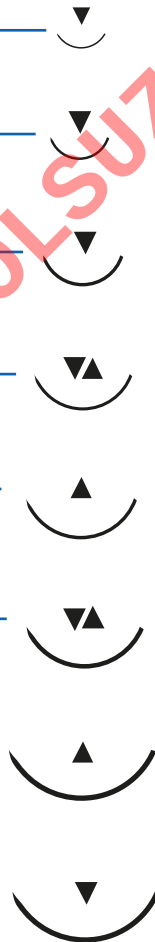


30 mm

35 mm

40 mm

50 mm


90 mm

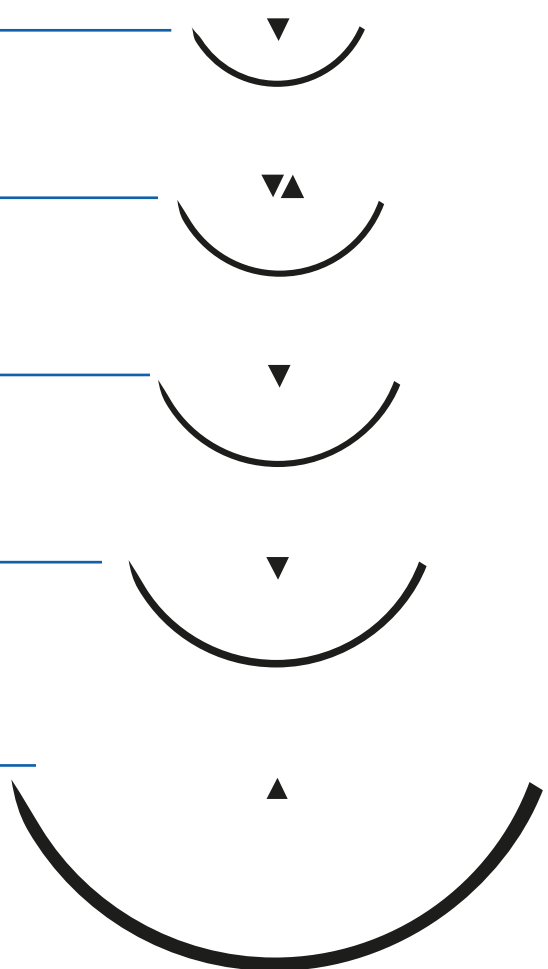


Special Needles

Loop


40 mm

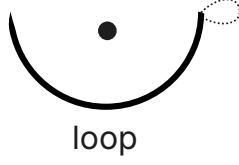

45 mm

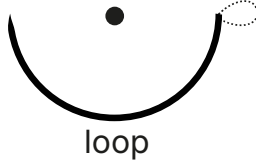

50 mm

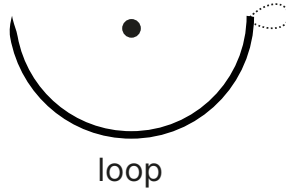


KONTROLSUZ KOPYADIR...



Absorbable Sutures

Absorbable Sutures

Wound Support*

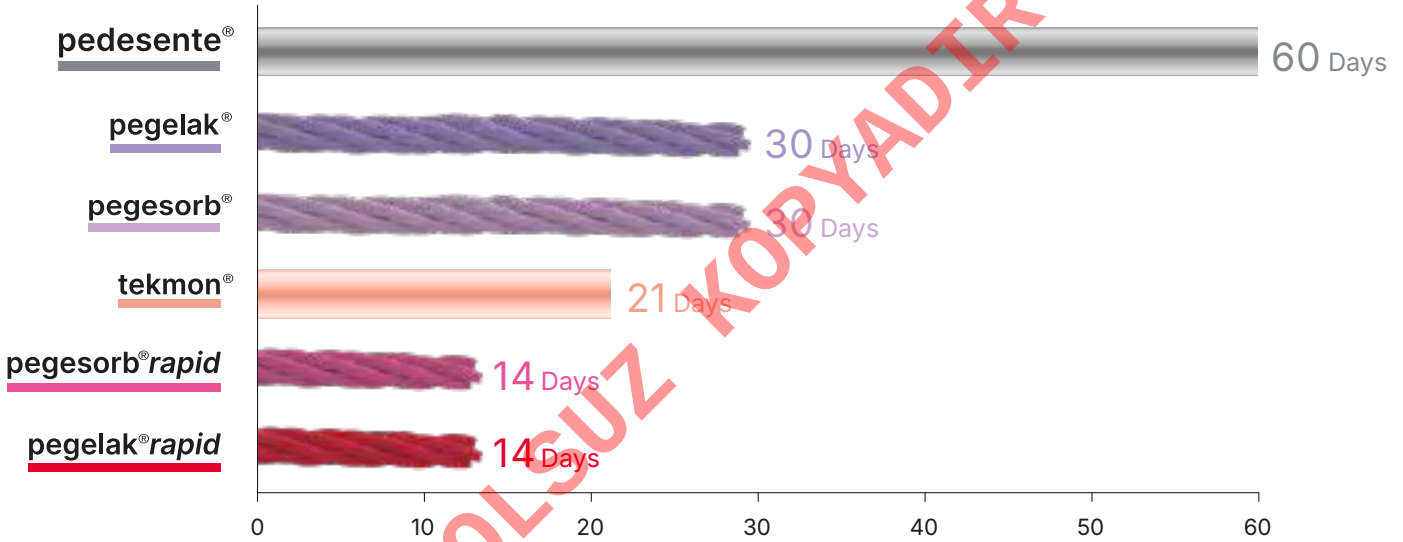

*Data is intended as guidance only.

Absorbable Sutures For All Your Needs



pegelak®





Description:	Absorbable, synthetic, braided, coated suture for medium term wound support with well proven record and composition
Characteristics:	Excellent handling properties Excellent knot security and tie-down Very high tensile strength Medium term wound support (~ 30 days)
Absorption:	Essentially complete in 56-70 days
Material :	polyglycolic acid (%90)-co-L-lactid acid (%10), PGLA [Copolymer of glycolide and L-lactide (90:10)]
Coating:	poly(glycolide-co-L-lactide) PGLA (30:70) and calcium stearate
Structure:	Braided (multifilament)
Colour :	Violet /undyed (beige)
Sterilization:	EO
Packaging:	Available in USP 7/0 (metric 0.5) through USP 2 (metric 5) in boxes of 12 (one dozen) aluminium pouches.

needle information		needle diameter	colour	suture length	0.5	0.7	1	1.5
					7/0	6/0	5/0	4/0
8 mm								
1/4 spatulated, double		280 μ	violet	45 cm		L1882		
		330 μ	violet	45 cm			L2882	
10 mm								
3/8 cutting		280 μ	violet	45 cm	L0042			
3/8 cutting, double		280 μ	violet	45 cm		L103		
13 mm								
3/8 round bodied		330 μ	undyed	45 cm		J131		
		330 μ	violet	45 cm		L131		
		330 μ	violet	60 cm			L231	
		330 μ	undyed	60 cm		J231		
		380 μ	violet	60 cm				L331
3/8 PL		280 μ	violet	45 cm		L139		
		330 μ	violet	45 cm			L239	
straight diamond		430 μ	violet	75 cm				L338

needle information		colour	suture length	1	1.5	2	3	3.5
				5/0	4/0	3/0	2/0	0
16 mm								
1/2 round bodied		violet	75 cm	L2160	L3160	L4160		
3/8 cutting		violet	75 cm	L2165	L3165	L4165		
		undyed	75 cm	J2165				
3/8 PL		violet	45 cm	L2169	L3169			
		undyed	75 cm	J2169	J3169			
18 mm								
1/2 round bodied		violet	75 cm		L3180	L4180		
straight diamond		violet	75 cm		L3188			
19 mm								
3/8 PL		violet	45 cm	L2199	L3199	L4199		
		undyed	45 cm	J2199	J3199	J4199		
20 mm								
1/2 round bodied		violet	75 cm		L3200	L4200		
3/8 cutting		violet	45 cm		L3203	L4203		
		violet	75 cm		L3213			
25 mm								
1/2 round bodied		violet	75 cm		L3250	L4250	L5250	
		violet	90 cm			L4260	L5240	L6260
1/2 cutting		violet	75 cm			L4254		
3/8 cutting		violet	45 cm		L3255	L4255	L5255	
3/8 PL		violet	45 cm		L3259	L4259		

needle information		colour	suture length	2	3	3.5	4	5
				3/0	2/0	0	1	2
26 mm								
1/2 round bodied		violet	75 cm	L4270	L5270			
3/8 cutting		violet	75 cm	L4263	L5263			
		undyed	75 cm	J4263	J5263			
30 mm								
1/2 round bodied		violet	75 cm	L4300	L5300	L6300	L7300	
		violet	90 cm		L5310	L6310	L7310	
3/8 cutting		violet	75 cm	L4303	L5303	L6303		
		undyed	75 cm	J4303	J5303	J6303		
35 mm								
1/2 round bodied		violet	75 cm		L5350	L6350	L7350	
		violet	90 cm				L7360	
3/8 cutting		violet	75 cm	L4353	L5353			
1/2 diamond		violet	90 cm		L5356	L6356		
36 mm								
5/8 round bodied		violet	75 cm			L6368		

needle information		colour	suture length	2	3	3.5	4	5
				3/0	2/0	0	1	2
40 mm								
1/2 round bodied		violet	75 cm		L5400	L6400	L7400	L8400
		violet	90 cm			L6410	L7410	
1/2 cutting		violet	90 cm			L6414	L7414	
1/2 diamond		violet	90 cm			L6406	L7406	

needle information		colour	suture length	1.5	2	3	3.5	4	5
				4/0	3/0	2/0	0	1	2

45 mm									
1/2 round bodied		violet	75 cm					L7450	L8450
		violet	90 cm					L7460	

50 mm									
1/2 round bodied		violet	75 cm			L5500	L6500	L7500	L8500
		violet	90 cm				L6510	L7510	L8510

60 mm									
straight cutting		violet	75 cm	L3608	L4608				
		undyed	75 cm	J3608					

needle information	colour	suture length	1.5	2	3	3.5	4	5
			4/0	3/0	2/0	0	1	2

63 mm								
3/8 blunt point liver needle		violet	75 cm					L7638
			90 cm					L7648

65 mm								
3/8 blunt point liver needle		violet	100 cm					L7658

76 mm								
1/2 round bodied, loop		violet	150 cm					L87603

	color	1	1.5	2	3	3.5	4	5
		5/0	4/0	3/0	2/0	0	1	2

without needle								
150 cm	violet	L25	L35	L45	L55	L65	L75	L85
6x45 cm	violet	L22	L32	L42	L52	L62	L72	L82
12x45 cm	violet	L24	L34	L44	L54	L64	L74	L84
17x45 cm	violet	L21	L31	L41	L51	L61	L71	L81
13x60 cm	violet	L23	L33	L43	L53	L63	L73	L83



pegelak®rapid





Description:	Absorbable, synthetic, braided, coated suture for short term wound support, ideal for fast healing wounds
Characteristics:	Quick absorption - increase in patient comfort and cosmetic results Excellent handling properties Excellent knot security and tie-down Short term wound support (~ 14 days)
Absorption:	Essentially complete in 42 days
Material :	polyglycolic acid (%90)-co-L-lactid acid (%10), PGLA [Copolymer of glycolide and L-lactide (90:10)]
Coating:	poly(glycolide-co-L-lactide) PGLA (30:70) and calcium stearate
Structure:	Braided (multifilament)
Colour :	Undyed (beige)
Sterilization:	Gamma Irradiation
Packaging:	Available in USP 6/0 (metric 0.7) through USP 2 (metric 5) in boxes of 12 (one dozen) aluminium pouches.

needle information		colour	suture length	0.7	1	1.5	2
				6/0	5/0	4/0	3/0
13 mm							
3/8 cutting		undyed	45 cm	Q133			
3/8 PL		undyed	45 cm	Q139	Q239		
straight diamond		undyed	75 cm			Q338	
16 mm							
1/2 round bodied		undyed	45 cm	Q1160			
		undyed	75 cm		Q2160	Q3160	Q4160
3/8 cutting		undyed	75 cm		Q2165	Q3165	Q4165
		undyed	75 cm			Q3163	
3/8 PL		undyed	75 cm		Q2169	Q3169	
18 mm							
1/2 round bodied		undyed	75 cm			Q3180	Q4180
19 mm							
3/8 PL		undyed	45 cm			Q3199	Q4199
20 mm							
1/2 round bodied		undyed	75 cm			Q3200	Q4200
3/8 cutting		undyed	75 cm			Q3213	Q4213
22 mm							
1/2 round bodied		undyed	45 cm			Q3220	Q4220

needle information	colour	suture length	1.5	2	3	3.5	4
			4/0	3/0	2/0	0	1

25 mm							
1/2 round bodied		undyed	75 cm	Q3250	Q4250	Q5250	
3/8 cutting		undyed	45 cm		Q4255	Q5255	
3/8 PL		undyed	45 cm	Q3259	Q4259		

26 mm							
3/8 cutting		undyed	75 cm		Q4263	Q5263	

30 mm							
1/2 round bodied		undyed	75 cm		Q4300	Q5300	Q6300 Q7300

35 mm							
1/2 round bodied		undyed	75 cm			Q5350	Q6350 Q7350
		undyed	90 cm				Q6360
3/8 cutting		undyed	75 cm		Q4353	Q5353	
1/2 diamond		undyed	90 cm			Q5356	Q6356

needle information	colour	suture length	3	3.5	4	5
			2/0	0	1	2


40 mm							
1/2 round bodied		undyed	75 cm	Q5400	Q6400	Q7400	Q8400
		undyed	90 cm	Q5410	Q6410	Q7410	
1/2 diamond		undyed	90 cm		Q6406	Q7406	

45 mm							
1/2 diamond		undyed	90 cm		Q6456		
1/2 round bodied		undyed	75 cm			Q7450	Q8450
		undyed	90 cm		Q6450		


50 mm							
1/2 round bodied		undyed	90 cm	Q5500	Q6500	Q7500	Q8500

	colour	1.5	2	3	3.5	4	5
		4/0	3/0	2/0	0	1	2

without needle							
150 cm	undyed	Q35	Q45	Q55	Q65	Q75	Q85
6x45 cm	undyed	Q32	Q42	Q52			

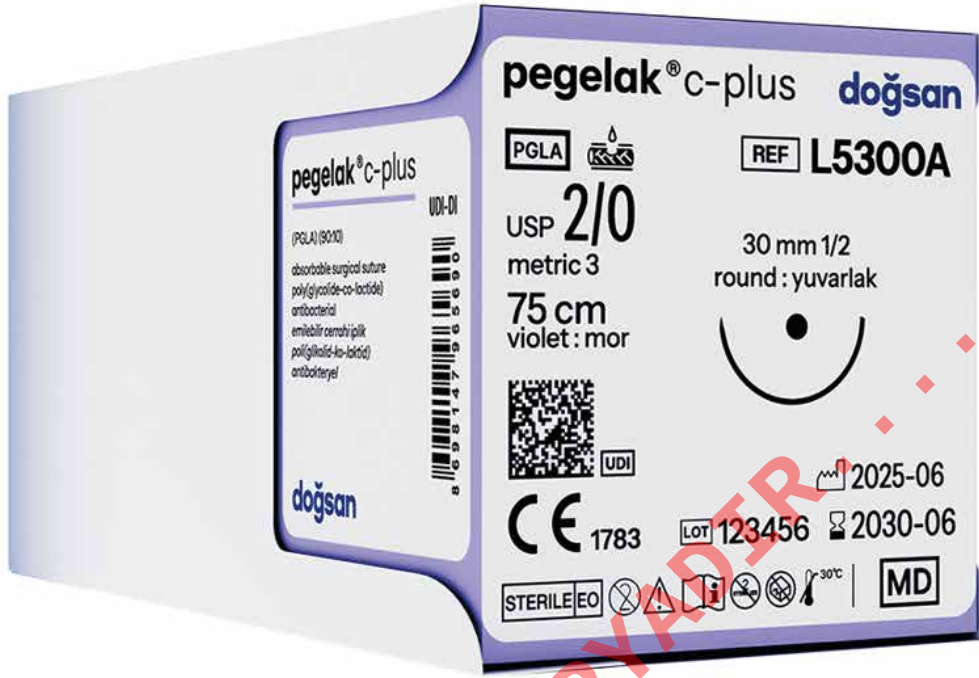

pegelak® C-plus

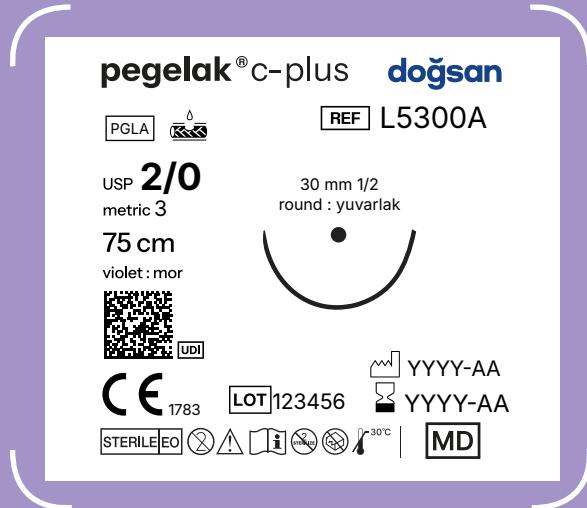

Pegelak[®] C-plus



Product Description:	Mid-term absorbable, sterile, synthetic, braided and coated antibacterial surgical suture made of PGLA (polyglactine 910 with Chlorhexidine Diacetate, a wide spectrum antibacterial. Antibacterial effect in addition to proven advantages of Pegelak suture (such as easy handling, high knot security, predictable absorption and high initial strength)
Antibacterial:	When used as intended, by reducing the possible surgical site infection formation along the suture against a wide spectrum of agents with chlorhexidine diacetate content.
Structure:	Multifilament, braided, coated
Polymer:	poly(glycolide-co-lactide (90:10) (PGLA/poliglactin 910) and chlorhexidine diacetate
Coating:	poly(glycolide-co-lactide) [PGLA (30:70)] and calcium stearate
Product range:	USP 6/0 (metric 0,7) – USP 2 (metric 5)
Color:	Violet or Undyed
Sterile Method:	Ethylene Oxide (EO)
Absorption:	Approximately 56-70 days. Wound Support: ≈ 30 days Tensile strength retention: First Day: 100 % After 2 weeks ~ 75% After 3 weeks ~ 40-50 %
Intended use:	Pegelak C-plus, is used soft tissue approximation and/or ligation. It is not used in cardiovascular Surgery, eye (ophthalmic) Surgery) and neurologic tissues. Pegelak C-plus antibacterial is indicated for use in general soft tissue approximation and/or ligation but not for use in cardiovascular, ophthalmic and neurological tissues. Being absorbable, it should not be used where extended approximation of tissue is required or should not be used in patients with known allergic reactions to chlorhexidine diacetate. Pegelak C-plus Sutures are recommended for a wide spectrum of applications such as, fascia, peritoneum, subcutaneous/intracutaneous tissue (fat), skin tissues and laparoscopic incisions.

Pegelak[®]C-plus

polyglycolic acid(%90)-co-lactic acid (%10), PGLA, antibacterial

needle information	needle diameter	color	suture length	0.7	1	1.5	2
				6/0	5/0	4/0	3/0

10 mm							
3/8 cutting		280 µ	violet	45 cm	L103A		

13 mm							
3/8 round bodied		330 µ	violet	45 cm	L131A		
				60 cm		L231A	
3/8 cutting		280 µ		45 cm	L133A		
		380 µ				L233A	

16 mm								
1/2 round bodied			violet	75 cm		L2160A	L3160A	L4160A
3/8 cutting						L2165A	L3165A	L4165A
3/8 cutting						L2163A	L3163A	L4163A
3/8 round bodied				45 cm		L2161A		

18 mm								
3/8 cutting			violet	75 cm			L3183A	L4183A

20 mm								
1/2 round bodied			violet	75 cm			L3200A	L4200A
3/8 cutting				45 cm			L3203A	L4203A
				75 cm			L3213A	L4213A

needle information	color	suture length	1.5	2	3	3.5	4
			4/0	3/0	2/0	0	1

25 mm							
1/2 round bodied		violet	75 cm	L3250A	L4250A	L5250A	
			90 cm		L4260A	L5240A	

26 mm							
3/8 cutting		violet	75 cm	L3263A	L4263A		

27 mm							
5/8 round bodied		violet	75 cm			L5278A	

30 mm							
1/2 round bodied		violet	75 cm		L4300A	L5300A	L6300A
			90 cm			L5310A	L7310A
1/2 cutting			75 cm			L5304A	
3/8 cutting					L4303A		L6303A

35 mm							
1/2 round bodied		violet	75 cm			L5350A	L6350A
			90 cm				L6360A

36 mm							
5/8 round bodied		violet	75 cm				L6368A

needle information		color	suture length	2	3	3.5	4	5
				3/0	2/0	0	1	2
40 mm								
1/2 round bodied		violet	75 cm		L5400A	L6400A	L7400A	L8400A
			90 cm			L6410A	L7410A	L8410A
1/2 cutting		violet	75 cm				L7404A	
			90 cm				L7414A	
45 mm								
1/2 round bodied		violet	75 cm				L7450A	L8450A
50 mm								
1/2 round bodied		violet	90 cm			L6510A	L7510A	L8510A
60 mm								
straight cutting		violet	75 cm	L4608A				
without needle								
without needle		violet	6x45 cm	L42A	L52A	L62A	L72A	



pegesorb®

doğsan

PGA



REF

G9672

USP 8/0

metric 0.4

45 cm

violet : mor

6 mm 3/8
spatulated - double :
spatül - çift



200µ



UDI

CE 1984

LOT

123456



2030-06

STERILE EO



30°C

MD

pegesorb®

UDI-DI

(PGA)
absorbable surgical suture
polyglycolic acid
emilabilir cerrahi iplik
poliglikolik asit



doğsan

pegesorb®

doğsan

pegesorb®

PGA



REF

G9672

USP 8/0

metric 0.4

45 cm

violet : mor

6 mm 3/8
spatuled - double:
spatül - çift



200 µ



UDI

CE 1984

LOT

123456



YYYY-AA



YYYY-AA

STERILE EO



30°C

MD



Description:	Absorbable, synthetic, braided, coated suture for medium term wound support with highest tensile strength
Characteristics:	Very high tensile strength Excellent handling properties Excellent knot security and tie-down Medium term wound support (~ 30 days)
Absorption:	Essentially complete in 60-90 days
Material :	polyglycolic acid, PGA %100
Coating:	poly(glycolide-co-L-lactide) PGLA (30:70) and calcium stearate (Also available with polycaprolactone & calcium stearate)
Structure:	Braided (multifilament)
Colour :	Violet /undyed (beige)
Sterilization:	EO
Packaging:	Available in USP 8/0 (metric 0.4) through USP 2 (metric 5) in boxes of 12 (one dozen) aluminium pouches

needle information		needle diameter	colour	suture length	0.4	0.5	0.7	1	1.5
					8/0	7/0	6/0	5/0	4/0
6 mm									
3/8 spatulated, double		200 μ	violet	45 cm	G9672				
8 mm									
1/4 spatulated, double		280 μ	violet	45 cm		G0882	G1882		
		330 μ	violet	45 cm				G2882	
3/8 cutting		200 μ	violet	45 cm		G083			
10 mm									
3/8 cutting		280 μ	violet	45 cm			G103		
3/8 cutting, double		280 μ	violet	45 cm			G1032		
13 mm									
3/8 round bodied		330 μ	undyed	45 cm			I131		
		330 μ	violet	45 cm			G131		
		330 μ	violet	60 cm				G231	
		380 μ	violet	60 cm					G331
3/8 PL		280 μ	violet	45 cm			G139		
		330 μ	violet	45 cm				G239	

needle information		colour	suture length	1	1.5	2	3
				5/0	4/0	3/0	2/0
16 mm							
1/2 round bodied		violet	75 cm	G2160	G3160	G4160	
3/8 round bodied		violet	45 cm	G2161			
3/8 cutting		violet	75 cm	G2165	G3165	G4165	
3/8 PL		violet	45 cm	G2169	G3169		
18 mm							
1/2 round bodied		violet	75 cm		G3180	G4180	
19 mm							
3/8 PL		violet	45 cm	G2199	G3199	G4199	
20 mm							
1/2 round bodied		violet	75 cm		G3200	G4200	
3/8 cutting		violet	45 cm		G3203	G4203	
		violet	75 cm		G3213		
		violet	90 cm			G4213	
25 mm							
1/2 round bodied		violet	75 cm		G3250	G4250	G5250
		violet	90 cm			G4260	G5260
3/8 cutting		violet	45 cm		G3255	G4255	G5255
3/8 PL		violet	45 cm		G3259	G4259	

needle information		colour	suture length	2	3	3.5	4	5
				3/0	2/0	0	1	2
26 mm								
3/8 cutting		undyed	75 cm	I4263	I5263			
		violet	75 cm	G4263	G5263			
1/2 round bodied		violet	75 cm	G4270	G5270			
27 mm								
5/8 round bodied		violet	75 cm		G5278			
30 mm								
1/2 round bodied		violet	75 cm	G4300	G5300	G6300	G7300	
1/2 cutting		violet	75 cm		G5304	G6304		
3/8 cutting		violet	75 cm	G4303	G5303	G6303		
35 mm								
1/2 round bodied		violet	75 cm	G4350	G5350	G6350	G7350	G8350
		violet	90 cm		G5360	G6360		
3/8 cutting		violet	75 cm	G4353	G5353			
1/2 diamond		violet	90 cm		G5356	G6356		
40 mm								
1/2 round bodied		violet	75 cm		G5400	G6400	G7400	G8400
			90 cm			G6410	G7410	G8410
1/2 cutting		violet	75 cm					G8404
			90 cm			G6414	G7414	G8414
1/2 diamond		violet	90 cm			G6406	G7406	G8406

needle information	colour	suture length	3	3.5	4	5
			2/0	0	1	2

45 mm

1/2 round bodied		violet	75 cm			G7450	G8450
		violet	90 cm			G7460	

50 mm

1/2 round bodied		violet	75 cm	G5500	G6500	G7500	G8500
		violet	90 cm		G6510	G7510	G8510

63 mm

3/8 blunt point liver needle	 blunt point	violet	75 cm			G7638	
------------------------------	--	--------	-------	--	--	-------	--

	colour	1	1.5	2	3	3.5	4	5
		5/0	4/0	3/0	2/0	0	1	2

without needle

150 cm	violet	G25	G35	G45	G55	G65	G75	G85
3x45 cm	violet	G28	G38	G48	G58	G68	G78	G88
6x45 cm	violet	G22	G32	G42	G52	G62	G72	G82
12x45 cm	violet	G24	G34	G44	G54	G64	G74	G84
17x45 cm	violet	G21	G31	G41	G51	G61	G71	G81



pegesorb® rapid





Description:	Absorbable, synthetic, braided, coated suture for short term wound support, ideal for fast healing wounds
Characteristics:	Quick absorption - increase in patient comfort, cosmetic results Excellent handling properties Excellent knot security and tie-down Short term wound support (~ 14 days)
Absorption:	Essentially complete in 42 days
Material :	polyglycolic acid, PGA 100%
Coating:	poly(glycolide-co-L-lactide) PGLA (30:70) and calcium stearate
Structure:	Braided (multifilament)
Colour :	Undyed (beige)
Sterilization:	Gamma Irradiation
Packaging:	Available in USP 5/0 (metric 1) through USP 2 (metric 5) in boxes of 12 (one dozen) aluminium pouches.

needle information		colour	suture length	1	1.5	2	3
				5/0	4/0	3/0	2/0
13 mm							
3/8 PL		undyed	45 cm	H239			
16 mm							
1/2 round bodied		undyed	75 cm	H2160	H3160	H4160	
3/8 cutting		undyed	75 cm	H2165	H3165	H4165	
3/8 PL		undyed	75 cm	H2169	H3169		
18 mm							
1/2 round bodied		undyed	75 cm		H3180	H4180	
19 mm							
3/8 PL		undyed	45 cm		H3199	H4199	
20 mm							
1/2 round bodied		undyed	75 cm		H3200	H4200	
3/8 cutting		undyed	75 cm		H3213	H4213	

needle information	colour	suture length	1.5	2	3	3.5	4
			4/0	3/0	2/0	0	1

25 mm							
1/2 round bodied		undyed	75 cm	H3250	H4250	H5250	
3/8 cutting		undyed	45 cm		H4255	H5255	
3/8 PL		undyed	45 cm	H3259	H4259		

26 mm							
3/8 cutting		undyed	75 cm		H4263	H5263	

30 mm							
1/2 round bodied		undyed	75 cm		H4300	H5300	H6300 H7300

35 mm							
1/2 round bodied		undyed	75 cm			H5350	H6350 H7350
			90 cm				H6360
3/8 cutting		undyed	75 cm		H4353	H5353	
1/2 diamond		undyed	90 cm			H5356	H6356

needle information	colour	suture length	3	3.5	4	5
			2/0	0	1	2

40 mm							
1/2 round bodied		undyed	75 cm	H5400	H6400	H7400	H8400
		undyed	90 cm	H5410	H6410	H7410	
1/2 diamond		undyed	90 cm		H6406	H7406	

45 mm							
1/2 round bodied		undyed	75 cm			H7450	H8450

50 mm							
1/2 round bodied		undyed	90 cm	H5500	H6500	H7500	H8500


	colour	1.5	2	3	3.5	4	5
		4/0	3/0	2/0	0	1	2

without needle							
150 cm	undyed	H35	H45	H55	H65	H75	H85



tekmon®



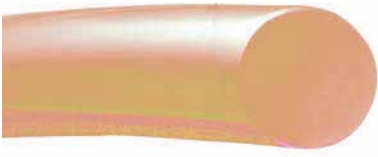

Description:	Absorbable, synthetic, monofilament suture for short-medium term wound support, with a smooth surface for fast healing wounds
Characteristics:	High initial strength Easy passage through tissue/smoothness of a monofilament Minimal tissue reaction Very good knot security Short -Medium term wound support (~ 21 days)
Absorption:	Essentially complete in 90-120 days
Material :	polyglycolide(%75)-co-caprolactone (%25), PGCL
Structure:	Monofilament
Colour :	Violet/Undyed (clear)
Sterilization:	EO
Packaging:	Available in USP 6/0 (metric 0.7) through USP 2 (metric 5) in boxes of 12 (one dozen) aluminium pouches.

needle information		needle diameter	colour	suture length	0.7	1	1.5	2
					6/0	5/0	4/0	3/0
8 mm								
1/4 spatulated, double		330 μ	undyed	45 cm	U1882	U2882		
10 mm								
3/8 cutting		330 μ	undyed	45 cm	U103			
13 mm								
3/8 PL		330 μ	undyed	45 cm	U139	U239		
3/8 round bodied		330 μ	undyed	60 cm		U231		
16 mm								
1/2 round bodied			undyed	75 cm		U2160		
			violet	75 cm			T3160	T4160
3/8 cutting			undyed	75 cm		U2165	U3165	U4165
			violet	90 cm			T3175	
3/8 PL			undyed	45 cm		U2169	U3169	
19 mm								
3/8 PL			undyed	45 cm		U2199	U3199	U4199

needle information	colour	suture length	1.5	2	3	3.5	4
			4/0	3/0	2/0	0	1

20 mm

1/2 round bodied		violet	75 cm	T3200	T4200			
3/8 cutting		undyed	75 cm	U3213	U4213			
		violet	45 cm	T3203	T4203			

25 mm


1/2 round bodied		violet	75 cm	T3250	T4250	T5250		
		undyed	75 cm		U4250	U5250		
3/8 cutting		violet	45 cm		T4255	T5255		
3/8 PL		undyed	45 cm	U3259	U4259			

26 mm

3/8 cutting		undyed	75 cm		U4263	U5263		
-------------	--	--------	-------	--	-------	-------	--	--


30 mm

1/2 round bodied		violet	75 cm		T4300	T5300	T6300	T7300
3/8 cutting		violet	75 cm			T5303	T6303	

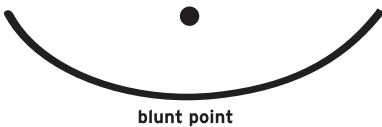needle information		colour	suture length	3	3.5	4	5
				2/0	0	1	2
35 mm							
1/2 round bodied		violet	75 cm	T5350	T6350	T7350	
3/8 cutting		violet	75 cm	T5353			
40 mm							
1/2 round bodied		violet	75 cm	T5400	T6400	T7400	T8400
		violet	90 cm		T6410	T7410	
1/2 cutting		violet	75 cm			T7404	T8404
45 mm							
1/2 round bodied		violet	75 cm			T7450	
50 mm							
1/2 round bodied		violet	75 cm	T5500		T7500	T8500
		violet	90 cm			T7510	
63 mm							
3/8 blunt point liver needle	 blunt point	violet	75 cm			T7638	
		colour	2	3	3.5	4	5
			3/0	2/0	0	1	2
without needle							
150 cm		violet	T45	T55	T65	T75	T85



pedesente®

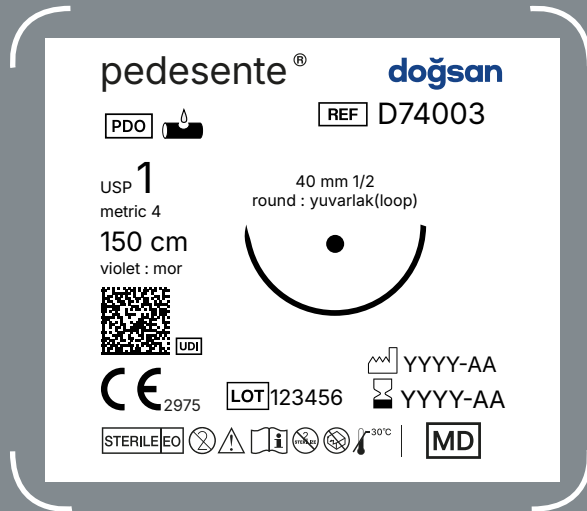




Description:	Absorbable, synthetic, monofilament suture for long term wound support, ideal for slow healing wounds
Characteristics:	Easy passage through tissue/smoothness of a monofilament High initial strength Slow decrease in tensile strength Minimal tissue reaction Ideal for long term support Long term wound support (~ 60 days)
Absorption:	Essentially complete in 180-210 days
Material :	polydioxanone, PDO 100%
Structure:	Monofilament
Colour :	Violet
Sterilization:	EO
Packaging:	Available in USP 7/0 (metric 0.5) through USP 2 (metric 5) in boxes of 12 (one dozen) aluminium pouches.


needle information		needle diameter	colour	suture lenght	0.5	0.7	1	1.5	2
					7/0	6/0	5/0	4/0	3/0
8 mm									
3/8 round bodied, double		220 μ	violet	75 cm	D0812				
		240 μ	violet	45 cm		D1812			
				75 cm		D1912			
3/8 KD, double		240 μ	violet	60 cm		D1872			
3/8 round bodied		240 μ	violet	30 cm		D181			
9 mm									
3/8 round bodied, double		200 μ	violet	75 cm	D0912				
10 mm									
3/8 cutting		280 μ	violet	45 cm		D103			
3/8 cutting, double		280 μ	violet	45 cm		D1032			
13 mm									
3/8 round bodied, double		220 μ	violet	45 cm	D0312				
		330 μ	violet	75 cm		D1312	D2312		
3/8 PL		280 μ	violet	45 cm		D139			
		330 μ	violet	45 cm			D239		
16 mm									
1/2 round bodied			violet	75 cm			D2160	D3160	D4160
3/8 round bodied			violet	75 cm			D2161	D3161	
1/2 round bodied, double			violet	90 cm			D21602	D31602	
3/8 cutting			violet	75 cm			D2165	D3165	D4165
3/8 PL			violet	45 cm			D2169	D3169	

needle information		colour	suture length	1	1.5	2	3	3.5
				5/0	4/0	3/0	2/0	0
19 mm								
3/8 PL		violet	45 cm	D2199	D3199	D4199		
20 mm								
1/2 round bodied		violet	75 cm		D3200	D4200		
3/8 cutting		violet	45 cm		D3203	D4203		
25 mm								
1/2 round bodied		violet	75 cm		D3250	D4250	D5250	D6250
3/8 cutting		violet	45 cm		D3255	D4255	D5255	
3/8 PL		violet	45 cm		D3259	D4259		
26 mm								
3/8 cutting		violet	75 cm			D4263		
30 mm								
1/2 round bodied		violet	75 cm			D4300	D5300	D6300
3/8 cutting		violet	75 cm				D5303	

needle information	colour	suture length	3	3.5	4	5
			2/0	0	1	2

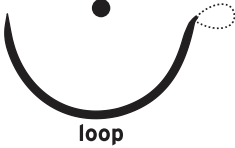40 mm							
1/2 round bodied		violet	75 cm	D5400	D6400	D7400	D8400
		violet	90 cm			D7410	D8410
1/2 round bodied, loop		violet	150 cm		D64003	D74003	

45 mm							
1/2 round bodied		violet	75 cm			D7450	D8450
1/2 round bodied, loop		violet	150 cm			D74503	

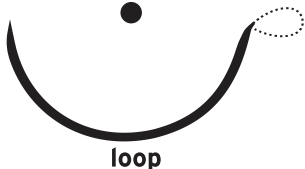50 mm							
1/2 round bodied		violet	75 cm		D6500	D7500	D8500
		violet	90 cm				D8510
1/2 round bodied, loop		violet	150 cm		D65003	D75003	

	propilen® polypropylene
	trofilen® polyvinylidene fluoride
	daylon® polyamide
	politer® polyester
	ipek/silk silk
	daxton-D polytetrafluoroethylene
	daxton-K expanded polytetrafluoroethylene
	tektel® 316L Stainless Steel Wire
	dopace® Twisted 316L stainless steel wire, with polyethylene (PE) coated (sheathed)

Non-Absorbable Sutures



propilen®





Description: Non-Absorbable, monofilament suture, made of specially extruded condensed polypropylene, very smooth and ideal for cardio vascular surgery

Characteristics: Ideal suture for cardio vascular surgery
Very smooth flow through the tissue
Very high tensile strength
High softness and flexibility
Minimum tissue reaction, ideal for infected tissues
Good knot security

Tensile strength: Permanent -Very high tensile strength

Material : polypropylene

Structure: Monofilament

Colour : Blue/Undyed (clear)

Sterilization: EO

Packaging: Available in USP 10/0 (metric 0.2) through USP 2 (metric 5) in boxes of 12 (one dozen) pouches.

needle information	needle diameter	suture length	0.2	0.3	0.4	0.5	0.7	1
			10/0	9/0	8/0	7/0	6/0	5/0

5 mm								
3/8 round bodied		100 μ	13 cm	M2205	M3205			

6 mm								
3/8 round bodied, double		160 μ	60 cm			P9612		
3/8 KD, double		160 μ	60 cm			P9672		

8 mm								
3/8 round bodied, double		160 μ	60 cm			P9812		
		220 μ	60 cm			P0812	P1812	
		180 μ	60 cm			P0822		
	 black	160 μ	60 cm			P9814		
		220 μ	60 cm			P0814	P1814	
		180 μ	60 cm			P9872		
3/8 KD, double		180 μ	60 cm			P9872		
		220 μ	60 cm			P0872	P1872	

9 mm								
3/8 round bodied, double		200 μ	75 cm			P0912		
	 black	200 μ	60 cm			P0914		

needle information		needle diameter	suture length	0.4	0.5	0.7	1
				8/0	7/0	6/0	5/0
10 mm							
3/8 round bodied, double		220 μ	60 cm		P0012		
		240 μ	60 cm			P1012	
			75 cm			P1112	
		330 μ	60 cm				P2012
	75 cm					P2112	
	 black	220 μ	60 cm		P0014		
		240 μ	75 cm			P1014	
	1/2 round bodied, double		220 μ	75 cm			P1002
330 μ			75 cm				P2002
3/8 KD, double		200 μ	60 cm		P0072		
		200 μ	60 cm		P0074		
	 black	280 μ	60 cm			P1074	
3/8 cutting		280 μ	45 cm			P103	

11 mm						
3/8 round bodied, double	 black	280 μ	60 cm			P1114

needle information		needle diameter	suture length	0.5	0.7	1	1.5	
				7/0	6/0	5/0	4/0	
13 mm								
3/8 round bodied		220 μ	60 cm	P031				
		280 μ	60 cm		P131			
			75 cm		P141			
		330 μ	60 cm			P231		
3/8 round bodied, double		220 μ	60 cm	P0312				
		280 μ	60 cm		P1312			
			75 cm					
		220 μ	75 cm	P0412				
		280 μ	75 cm		P1412			
			75 cm			P2312		
		380 μ	90 cm			P2412		
		380 μ	75 cm				P3312	
			220 μ	60 cm	P0314			
			240 μ	60 cm		P1314		
75 cm				P1414	P2314			
1/2 round bodied		280 μ			P130			
		330 μ				P230		
1/2 round bodied, double		280 μ	75 cm		P1302			
		330 μ	75 cm			P2302		
3/8 KD, double		280 μ	75 cm		P1372			
		330 μ	75 cm			P2372		
3/8 cutting		280 μ	45 cm		P133			
		380 μ	45 cm			P233	P333	

needle information	colour	suture length	0.7	1	1.5	2
			6/0	5/0	4/0	3/0

16 mm							
1/2 round bodied		blue	60 cm	P1160			
		blue	75 cm		P2160	P3160	P4160
1/2 round bodied, double		blue	90 cm		P21602	P31602	P41602
		black	90 cm		P21604	P31604	
1/2 KD, double		blue	90 cm		P21672	P31672	P41672
3/8 round bodied, double		blue	90 cm		P21612	P31612	
3/8 cutting		blue	45 cm	P1165	P2165	P3165	
		blue	75 cm				P4165
3/8 cutting		blue	45 cm	P1163	P2163	P3163	
		blue	75 cm				P4163
		blue	75 cm		P2173	P3173	
3/8 PL		blue	45 cm	P1169		P3169	

18 mm							
1/2 round bodied		blue	75 cm			P3180	P4180
		undyed	75 cm			X3180	
3/8 cutting		blue	75 cm		P2183	P3183	
		undyed	75 cm			X3183	

19 mm							
3/8 PL		blue	45 cm		P2199	P3199	P4199

needle information	colour	suture length	1	1.5	2	3	3.5
			5/0	4/0	3/0	2/0	0

20 mm								
1/2 round bodied		blue	75 cm	P2200	P3200	P4200	P5200	
1/2 round bodied, double		blue	90 cm	P2200 2	P3200 2	P4200 2		
	 black	blue	90 cm		P3200 4			
1/2 KD, double		blue	90 cm	P2207 2	P3207 2			
3/8 cutting		blue	45 cm	P2205	P3205	P4205		
3/8 cutting		blue	45 cm	P2203	P3203	P4203		
		blue	75 cm	P2213	P3213	P4213		

22 mm								
1/2 round bodied		blue	90 cm	P2220 2	P3220 2			
1/2 cutting		blue	45 cm			P4224		

25 mm								
1/2 round bodied		blue	75 cm		P3250	P4250	P5250	P6250
		undyed	75 cm			X4250		
1/2 round bodied, double		blue	90 cm		P3250 2	P4250 2	P5250 2	
	 black	blue	90 cm		P3250 4	P4250 4		
3/8 PL		blue	45 cm	P2259	P3259	P4259		
1/2 KD, double		blue	90 cm		P3257 2	P4257 2	P5257 2	
3/8 cutting		blue	45 cm		P3255	P4255	P5255	
		blue	75 cm		P3245			

needle information	colour	suture length	2	3	3.5	4
			3/0	2/0	0	1

26 mm						
3/8 cutting		blue	75 cm	P4263	P5263	
		undyed	75 cm	X4263		

30 mm						
1/2 round bodied		blue	75 cm	P4300	P5300	P6300
		blue	100 cm			P7300
1/2 round bodied, double		blue	90 cm	P43002	P53002	
1/2 cutting		blue	75 cm		P5304	
3/8 cutting		blue	75 cm	P4303	P5303	P6303


35 mm						
1/2 round bodied		blue	75 cm		P5350	
		blue	100 cm			P6350 P7350
3/8 cutting		blue	75 cm	P4353	P5353	

needle information	suture length	3	3.5	4	5
		2/0	0	1	2


40 mm

1/2 round bodied		100 cm	P5400	P6400	P7400	P8400
1/2 cutting		100 cm		P6404	P7404	
3/8 cutting		100 cm	P5403		P7403	
3/8 cutting		100 cm	P5405			

45 mm

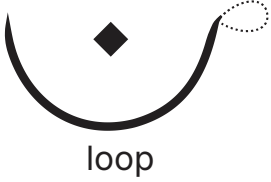1/2 round bodied	 loop	150 cm		P64503	P74503	
1/2 diamond		100 cm			P7456	
1/2 diamond	 loop	200 cm			P74563	

needle information	suture length	3	3.5	4	5
		2/0	0	1	2

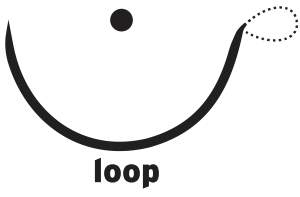50 mm					
1/2 round bodied		100 cm			P7500
	 loop	150 cm		P6500 3	P75003
		200 cm			P75103
1/2 cutting		100 cm			P7504

needle information	suture length	2	3	3.5	4
		3/0	2/0	0	1

60 mm					
straight cutting		75 cm	P4608	P5608	



trofilen® **doğsan**

PVDF  **REF F133**

USP **6/0**
metric 0.7

45 cm
blue : mavi

13 mm 3/8
cutting: keskin


280µ

 **UDI**

CE 1984 **LOT** 123456  2025-06  2030-06

STERILE **EO**  **MD**

trofilen®



trofilen® **doğsan**

PVDF  **REF F133**

USP **6/0**
metric 0.7

45 cm
blue : mavi

13 mm 3/8
cutting : keskin


280 µ

 **UDI**

CE 1984 **LOT** 123456  YYYY-AA  YYYY-AA

STERILE **EO**  **MD**



- Description:** Non-Absorbable, monofilament suture, made of specially extruded Polyvinylidene fluoride (PVDF) a very inert material ideal for plastic and cardiovascular surgery
- Characteristics:** Very smooth flow through the tissue
Highest tensile strength in monofilaments
Excellent pliability and flexibility
Minimum tissue reaction, ideal in infected tissues
- Tensile strength:** Permanent - Very high tensile strength
- Material :** polyvinylidene fluoride (PVDF)
- Structure:** Monofilament
- Colour :** Blue
- Sterilization:** EO
- Packaging:** Available in USP 8/0 (metric 0.4) through USP 2 (metric 5) in boxes of 12 (one dozen) pouches.

needle information	needle diameter	suture length	0.4	0.5	0.7	1
			8/0	7/0	6/0	5/0

8 mm							
3/8 round bodied, double		160 μ	60 cm	F9812			
		220 μ	60 cm		F0812	F1812	
	 black	160 μ	60 cm	F9814			
		220 μ	75 cm		F0814		
		220 μ	60 cm			F1814	
3/8 KD, double		180 μ	60 cm	F9872			
		220 μ	60 cm		F0872	F1872	
3/8 cutting		220 μ	45 cm		F083		

10 mm							
3/8 round bodied, double		220 μ	60 cm		F0012		
		240 μ	60 cm			F1012	
3/8 KD, double	 black	220 μ	60 cm		F0014		
		220 μ	60 cm		F0072		
		330 μ	60 cm				F2072
	 black	200 μ	60 cm		F0074		
3/8 cutting		280 μ	45 cm			F103	

needle information		needle diameter	suture length	0.5	0.7	1	1.5
				7/0	6/0	5/0	4/0
13 mm							
3/8 round bodied		280 μ	60 cm		F131		
		330 μ	60 cm			F231	
3/8 round bodied, double		220 μ	60 cm	F0312			
		280 μ	60 cm		F1312		
			75 cm		F1412		
	330 μ	75 cm			F2312		
	 black	220 μ	60 cm	F0314			
		240 μ	60 cm		F1314		
		240 μ	75 cm		F1414		
330 μ		75 cm			F2314		
1/2 round bodied, double		330 μ	75 cm			F2302	
3/8 cutting		280 μ	45 cm		F133		
		380 μ	45 cm			F233	F333

needle information		suture length	0.7	1	1.5	2
			6/0	5/0	4/0	3/0
16 mm						
1/2 round bodied		60 cm	F1160			
		75 cm		F2160	F3160	F4160
1/2 round bodied, double	 black	90 cm		F21602	F31602	F41602
		90 cm		F21604	F31604	
1/2 KD, double		90 cm		F21672	F31672	F41672
3/8 cutting		45 cm	F1165	F2165	F3165	
		75 cm				F4165
3/8 cutting		45 cm	F1163	F2163	F3163	
		75 cm				F4163
3/8 PL		45 cm	F1169	F2169	F3169	
19 mm						
3/8 PL		45 cm		F2199	F3199	F4199

needle information	suture length	1	1.5	2	3	3.5
		5/0	4/0	3/0	2/0	0

20 mm							
1/2 round bodied		75 cm	F2200	F3200	F4200	F5200	
1/2 round bodied, double		90 cm	F22002	F32002	F42002		
	 black	90 cm		F32004			
1/2 KD, double		90 cm	F22072	F32072			
3/8 cutting		45 cm	F2205	F3205	F4205		
3/8 cutting		45 cm	F2203	F3203	F4203		


22 mm							
1/2 round bodied		75 cm	F2220	F3220	F4220	F5220	


25 mm							
1/2 round bodied		75 cm		F3250	F4250	F5250	
1/2 round bodied, double		90 cm		F32502	F42502	F52502	
	 black	90 cm		F32504	F42504		
3/8 PL		45 cm	F2259	F3259	F4259		
1/2 KD, double		90 cm		F32572	F42572	F52572	
3/8 cutting		45 cm		F3255	F4255	F5255	

26 mm							
3/8 cutting		75 cm			F4263	F5263	

needle information		suture length	2	3	3.5	4	5
			3/0	2/0	0	1	2
30 mm							
1/2 round bodied		75 cm	F4300	F5300	F6300		
		100 cm				F7300	
1/2 cutting		75 cm		F5304			
3/8 cutting		75 cm	F4303	F5303	F6303		
35 mm							
1/2 round bodied		75 cm		F5350			
		100 cm			F6350	F7350	
3/8 cutting		75 cm	F4353	F5353			
3/8 cutting		75 cm				F7355	
36 mm							
5/8 round bodied		75 cm				F7368	
40 mm							
1/2 round bodied		100 cm		F5400	F6400	F7400	
1/2 cutting		100 cm			F6404	F7404	
3/8 cutting		100 cm		F5403		F7403	

needle information	suture length	2	3	3.5	4	5
		3/0	2/0	0	1	2

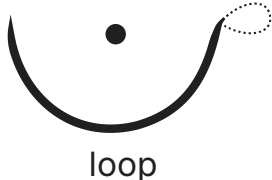45 mm						
1/2 round bodied	 loop	150 cm			F64503	F74503

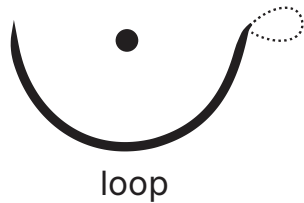50 mm						
1/2 round bodied	  loop	100 cm				F7500
		150 cm			F65003	F75003
1/2 cutting		100 cm				F7504


60 mm						
straight cutting		75 cm	F4608	F5608		

90 mm						
3/8 cutting, retention suture		2x50 cm				F79056



daylon®



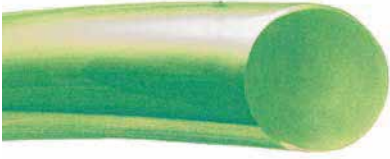

Description:	Non-Absorbable, monofilament suture, made of polyamide classic suture for plastic and ophthalmic surgery
Characteristics:	Uniform diameter Minimal memory High Flexibility Very good knot security Smooth passage through tissue
Tensile strength:	Permanent
Material :	Polyamide 6 -6.6
Structure:	Monofilament
Colour :	Blue : USP 6/0 - USP 2 Black: USP 11 /0 - USP 8/0
Sterilization:	EO
Packaging:	Available in USP 11 /0 (metric 0.01) through USP 2 (metric 5) in boxes of 12 (one dozen) pouches.

needle information	needle diameter	colour	suture length	0.1	0.2	0.3	0.4
				11/0	10/0	9/0	8/0



4 mm							
3/8 round bodied		50 µ	black	13cm			
		70 µ	black	13cm	M1102		
		100 µ	black	13cm		M2102	M3102
		140 µ	black	13cm			M4112


5 mm							
3/8 round bodied		100 µ	black	13cm		M2105	M3105
		140 µ	black	13cm			M4115
3/8 spatulated, double		200 µ	black	30 cm		M2185	

6 mm							
3/8 spatulated, double		150 µ	black	30 cm		M2177	M3167
		200 µ	black	45 cm		M2187	
3/8 round bodied		140 µ	black	13cm			M4117

needle information		needle diameter	colour	suture length	0.7	1	1.5	2
					6/0	5/0	4/0	3/0
10 mm								
3/8 cutting		280 μ	blue	45 cm	Y103			
13 mm								
3/8 cutting		280 μ	blue	45 cm	Y133			
3/8 PL		280 μ	blue	45 cm	Y139			
		330 μ	blue	45 cm		Y239		
16 mm								
3/8 cutting		blue	45 cm	Y1163	Y2163	Y3163		
		blue	75 cm					Y4163
		undyed	75 cm					V4163
3/8 PL		blue	45 cm	Y1169	Y2169	Y3169		
19 mm								
3/8 PL		blue	45cm		Y2199	Y3199	Y4199	
		undyed	45cm					V4199
20 mm								
3/8 cutting		blue	45cm		Y2203	Y3203	Y4203	
3/8 cutting		blue	45cm			Y3205	Y4205	

needle information		colour	suture length	1	1.5	2	3	3.5	4
				5/0	4/0	3/0	2/0	0	1
25 mm									
3/8 cutting		blue	45 cm		Y3255	Y4255			
3/8 PL		blue	45 cm	Y2259	Y3259	Y4259			
26 mm									
3/8 cutting		blue	75 cm			Y4263	Y5263		
		undyed	75 cm			V4263	V5263		
30 mm									
1/2 round bodied		blue	100 cm					Y6300	Y7300
3/8 cutting		blue	75 cm			Y4303	Y5303		
		undyed	75 cm				V5303		
35 mm									
1/2 round bodied		blue	100 cm					Y6350	Y7350
40 mm									
1/2 round bodied		blue	100 cm				Y5400	Y6400	Y7400
1/2 cutting		blue	100 cm					Y6404	Y7404
3/8 cutting		blue	100 cm						Y7403

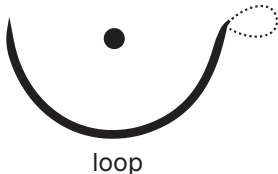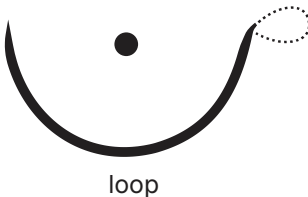needle information	colour	suture length	2	3	3.5	4
			3/0	2/0	0	1
45 mm						
1/2 round bodied	 loop	blue	150 cm		Y64503	Y74503
50 mm						
1/2 round bodied	 loop	blue	150 cm		Y65003	Y75003
60 mm						
straight cutting		blue	75 cm	Y4608	Y5608	

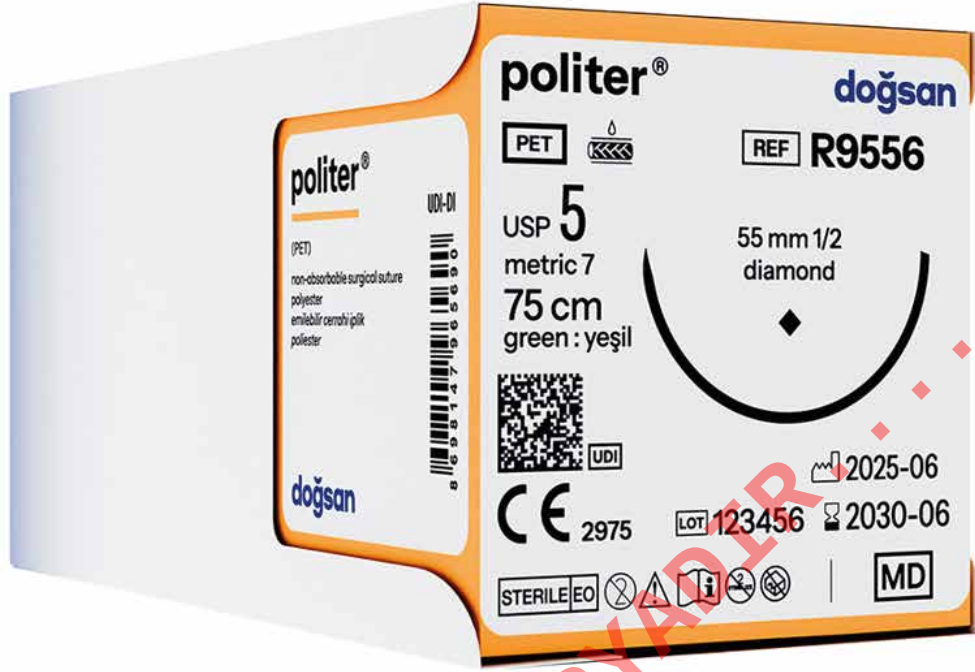

politer®





Description:	Non-Absorbable, synthetic, braided, coated suture, made by braiding of fine polyester filaments, ideal for cardiovascular & implant surgery for very long term support
Characteristics:	Very high tensile strength Permanent support Excellent knot security Excellent handling Providing maximum security under stress
Tensile strength:	permanent, no known significant change in tensile strength retention in vivo.
Material:	polyethylene terephthalate (polyester)
Coating:	silicone based
Structure:	Braided (multifilament)
Colour:	Green / White
Sterilization:	EO
Packaging:	Available in USP 6/0 (metric 1) through USP 5 (metric 7) in boxes of 12 (one dozen) pouches. Or in the case of multistrand sutures with PTFE Pledgets (K-PAK), half a dozen (6) boxes.

needle information		needle diameter	colour	suture length	1	1.5	2	3
					5/0	4/0	3/0	2/0
8 mm								
1/4 spatulated, double		330 μ	white	45 cm	O2882			
		330 μ	green	45 cm	R2882			
13 mm								
3/8 cutting		380 μ	green	45 cm	R233	R333		
3/8 round bodied, double		330 μ	green	75 cm	R2312			
		380 μ	green	75 cm		R3312		
16 mm								
1/2 round bodied			green	75 cm		R3160	R4160	
1/2 round bodied, double			green	90 cm		R31602	R41602	
3/8 cutting			green	75 cm		R3165	R4165	
17 mm								
1/2 round bodied, double			green	90 cm		R31702	R41702	R51702
1/2 round bodied, double	  With PTFE Pledget 6 x 3 mm		green	75 cm				R5170P
			white	75 cm				O5170P
1/2 diamond, double			green	90 cm	R21762	R31762	R41762	R51762
1/2 diamond, double	  With PTFE Pledget 6 x 3 mm		green	75 cm				R5176P
18 mm								
3/8 diamond, double			green	75 cm			R41862	

needle information	needle diameter	colour	suture length	1	1.5	2	3
				5/0	4/0	3/0	2/0

20 mm								
1/2 round bodied, double			green	90 cm		R3200 2	R42002	R52002
			green	75 cm		R31902	R41902	R51902
1/2 round bodied, double			green	75 cm				R5190P
			white	75 cm				05190P
	 With PTFE Pledget 6 x 3 mm							

needle information	colour	suture length	2	3	3.5	4	5
			3/0	2/0	0	1	2

25 mm								
1/2 round bodied			green	75 cm	R4250	R5250	R6250	
1/2 round bodied, double			green	75 cm	R4240 2	R52402		
			green	90 cm	R4250 2	R52502	R62502	
	  With PTFE Pledget 6 x 3 mm		green	75 cm		R5240P		
			white	75 cm		O5240P		
1/2 diamond, double			green	90cm	R42562	R52562		
			green	75cm		R5256P		
	  With PTFE Pledget 6 x 3 mm		white	75cm		O5256P		
3/8 cutting			green	75 cm	R4255	R5255		

30 mm								
1/2 round bodied		green	75 cm		R5300	R6300		
3/8 cutting		green	75 cm	R4303	R5303			

needle information	colour	suture length	2	3	3.5	4	5
			3/0	2/0	0	1	2

35 mm							
1/2 round bodied		green	75 cm			R6350	R7350
3/8 cutting		green	75 cm				R7355

37 mm							
1/2 round cutting		green	75 cm				R7374

40 mm								
1/2 round bodied		green	75 cm			R6400	R7400	R8400

45 mm								
1/2 diamond		green	75 cm					R8456

needle information	suture length	3	3.5	4	5	7
		2/0	0	1	2	5

50 mm							
1/2 round bodied		75 cm			R7500	R8500	
straight cutting		75 cm	R5508				

needle information	suture length	3	3.5	4	5	7
		2/0	0	1	2	5

55 mm						
1/2 diamond		75 cm				R9556

90 mm						
3/8 cutting		100 cm			R7905	

	2	3	3.5	4
	3/0	2/0	0	1

without needle				
17x45 cm	R41	R51	R61	R71
13x60 cm	R43	R53	R63	



K3-PAK®

needle information		colour	suture length	3	3.5	4	5
				2/0	0	1	2
17 mm							
1/2 round bodied	 With PTFE Pledget 6 x 3 mm	5 X green 5 X white	10x75 cm	R5160K			
25 mm							
1/2 round bodied	 With PTFE Pledget 6 x 3 mm	5 X green 5 X white	10x75 cm	R5240K			
		5 X green 5 X white	10x90 cm	R5250K			



ipek/silk





- Description:** Non-Absorbable, braided, coated suture, made by braiding of fine processed natural silk fibers called fibroins, known for handling properties and used in many branches of surgery
- Characteristics:** Excellent handling
Excellent knot security
Coated suture for smooth flow through the tissue
- Tensile strength:** progressive degradation of the silk fibers resulting in gradual loss of tensile strength
- Material :** natural silk fibers of Bombyx mori
- Coating:** silicone based
- Structure:** Braided (multifilament) /twisted virgin silk for USP 8/0
- Colour :** Black / Blue only for 8/0 virgin silk
- Sterilization:** EO
- Packaging:** Available in USP 8/0 (metric 0.4) through USP 3 (metric 6) in boxes of 12 (one dozen) pouches.

needle information	needle diameter	suture length	0.4	0.5	0.7
			8/0	7/0	6/0
6 mm					
3/8 spatulated, double		140 μ	30 cm	M4367	
		200 μ	45 cm	M4387	

needle information		needle diameter	suture length	0.7 	1	1.5
				6/0	5/0	4/0
8 mm						
1/4 spatulated, double		280 μ	45 cm	S1882		
		330 μ	45 cm		S2882	
10 mm						
3/8 round bodied		330 μ	75 cm		S201	
3/8 cutting		330 μ	45 cm	S103		
13 mm						
3/8 round bodied		330 μ	45 cm	S131		
		330 μ	75 cm		S231	
		380 μ	75 cm			S331
3/8 cutting		280 μ	45 cm	S133		
		380 μ	45 cm		S233	S333

needle information	suture length	0.7	1	1.5	2	3
		6/0	5/0	4/0	3/0	2/0

16 mm						
1/2 round bodied		45 cm	S1160			
		75 cm		S2160	S3160	S4160
3/8 cutting		45 cm	S1165			
		75 cm		S2165	S3165	◆ S4165

18 mm						
1/2 round bodied		75 cm		S2180	S3180	S4180
3/8 cutting		75 cm		S2185	S3185	S4185

20 mm							
1/2 round bodied		75 cm		S2200	S3200	S4200	S5200
3/8 cutting		75 cm		S2205	S3205	S4205	
3/8 cutting		75 cm		S2203	S3203	S4203	S5203

22 mm							
1/2 cutting		75 cm			S3224	S4224	S5224

25 mm							
3/8 cutting		75 cm		S2255	S3255	S4255	S5255

needle information	suture length	1.5	2	3	3.5	4
		4/0	3/0	2/0	0	1

26 mm

1/2 round bodied		75 cm	S3260	S4260	S5260	S6260	
3/8 cutting		75 cm	S3263	S4263	S5263	S6263	

30 mm

1/2 round bodied		75 cm		S4300	S5300	S6300	S7300
3/8 cutting		75 cm		S4303	S5303	S6303	
1/2 cutting		75 cm			S5304		S7304

35 mm

1/2 round bodied		75 cm			S5350	S6350	S7350
3/8 cutting		75 cm					S7355
3/8 cutting		75 cm		S4353	S5353		

37 mm

1/2 cutting		75 cm				S6374	S7374
1/2 round bodied		75 cm				S6370	S7370

needle information		suture length	2	3	3.5	4	5	6
			3/0	2/0	0	1	2	3
40 mm								
1/2 round bodied		75 cm		S5400	S6400	S7400	S8400	
3/8 cutting		75 cm			S6403	S7403	S8403	
1/2 cutting		75 cm				S7404	S8404	
50 mm								
straight cutting		75 cm	S4508	S5508				
3/8 cutting		75 cm						S9503
1/2 cutting		75 cm				S7500		S9500
60 mm								
straight cutting		75 cm	S4608	S5608				
without needle		1.5	2	3	3.5	4	5	6
		4/0	3/0	2/0	0	1	2	3
17x45 cm	S31	S41	S51	S61	S71	S81		
13x60 cm	S33	S43	S53	S63	S73	S83		
180 cm	S32	S42	S52	S62	S72	S82	S92	
10x75 cm		S46	S56	S66	S76	S86		



daxton®-D





Description:	Daxton D is a non-absorbable, sterile surgical monofilament suture, composed of polytetrafluoroethylene (PTFE). PTFE suture has nonporous structure.
Characteristics:	Compatible needle/suture diameters, excellent tensile strength, secure knot positioning, minimum plaque formation, minimum tissue reaction, easy handling, high softness and minimal surface friction, elasticity.
Tensile strength:	Permanent
Material:	Polytetrafluoroethylene (PTFE)
Structure:	Monofilament
Colour:	White
Sterilization:	EO
Packaging:	Available in USP7/0 (metric 0.5) through USP 2/0 (metric 3) in boxes of 12 (one dozen) pouches.

needle information		colour	suture length	1	1.5	2	3
				5/0	4/0	3/0	2/0
13 mm							
1/2 round bodied		white	75 cm		XD330		
3/8 cutting		white	75 cm	XD233	XD333		
16 mm							
3/8 cutting		white	75 cm	XD2163	XD3163	XD4163	
3/8 PL		white	75 cm			XD4169	
18 mm							
3/8 cutting		white	75 cm	XD2183	XD3183	XD4183	
19 mm							
3/8 cutting		white	75 cm			XD4193	XD5193
50 mm							
straight round		white	90 cm				XD5509
60 mm							
straight cutting		white	75 cm				XD5608



daxton®-K







Description:	daxton®-K is a non-absorbable, sterile surgical monofilament suture, composed of polytetrafluoroethylene (PTFE), with chemical formula (C ₂ F ₄) _n , that has been expanded to an ePTFE suture with a porous microstructure and contains air.
Characteristics:	Compatible needle/suture diameters, excellent tensile strength, secure knot positioning, minimum tissue reaction, easy handling, high softness and minimal surface friction, elasticity.
Tensile strength:	Permanent
Material:	expanded Polytetrafluoroethylene (ePTFE)
Structure:	Monofilament
Colour:	White
Sterilization:	EO
Packaging:	Available in KD01 through KD5 in boxes of 12 (one dozen) pouches.

needle information		colour	suture length	KD-2	KD-3	KD-4	KD-5
13 mm							
1/2 round bodied		white	75 cm				XK2302
16 mm							
1/2 round bodied		white	75 cm			XK31602	
			90 cm			XK21602	
18 mm							
1/2 round bodied		white	75 cm		XK41802		
						XK31802	
3/8 cutting			45 cm				XK2183
22 mm							
1/2 round bodied		white	90 cm			XK32202	
25 mm							
1/2 round bodied		white	90 cm		XK42502		
						XK32502	
35 mm							
1/2 round bodied		white	90 cm	XK53502			



tektel®

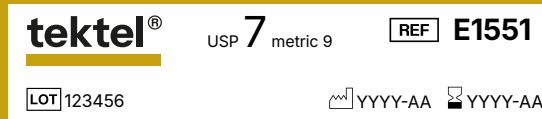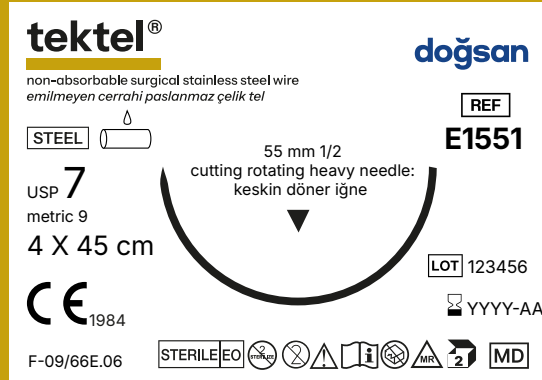





Description:	Non-absorbable stainless steel wire
Characteristics:	360 ° rotating or fixed head with needle alternatives (Cutting, diamond, trocar point), very high tensile strength, minimum tissue reaction, easy pliability and handling.
Tensile strength:	Permanent
Material:	Monofilament
Structure:	316L
Colour:	Metallic grey
Sterilization:	EO
Packing:	Available in USP 2 (metric 5) - USP 8 (metric 10) in boxes of 2 pouches with rotating or fixed head needle alternatives

needle information	suture length	5	6	7	8	9
		2	3&4	5	6	7


45 mm						
1/2 diamond	 fixed head needle	4 X 45 cm	E7457	E8457		

48 mm						
1/2 diamond	 rotating needle ↻	4 X 45 cm			E9486	E0486
		2 X 75 cm			E9488	
	 fixed head needle	4 X 45 cm				E0487

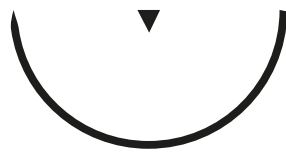55 mm						
1/2 cutting	 rotating needle ↻	2 X 75 cm			E9553	E1553
		4 X 45 cm			E9551	E0551
	 fixed head needle	4 X 45 cm				E0554

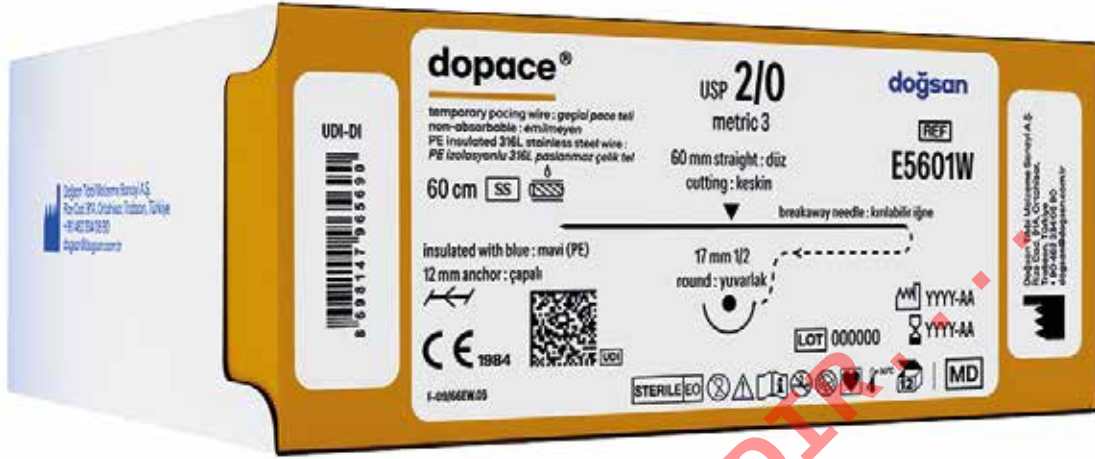




dopace®

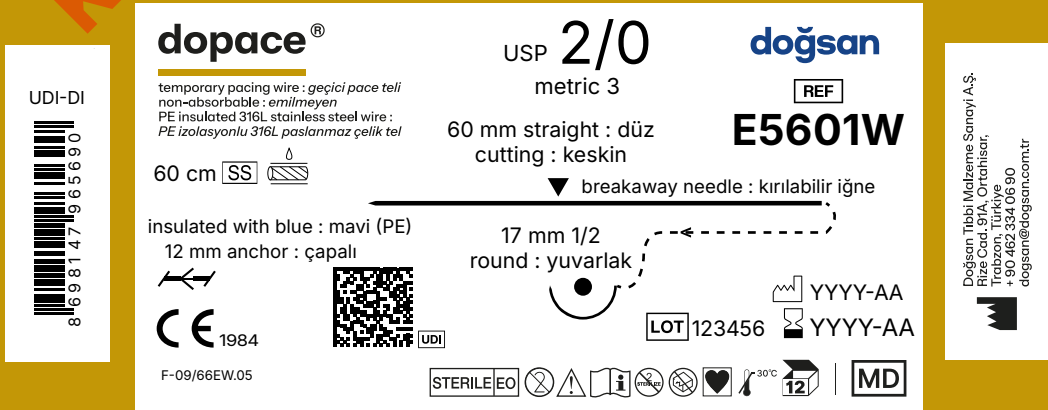



Description:	Temporary pacing wire, PE insulated 316L stainless steel wire
Characteristics:	Intracorporeal curved needle and extracorporeal straight breakaway needle, with or without anchor. With anchor option no need for suture to fix it and less trauma on myocardium muscle.
Tensile strength:	Permanent
Material:	Twisted 316L stainless steel wire, with polyethylene (PE) coated (sheathed)
Structure:	-
Colour:	-
Sterilization:	EO
Packing:	Available in USP 2/0 (metric 3) in boxes of 12 (one dozen) pouches.

needle information	suture length	3
		2/0

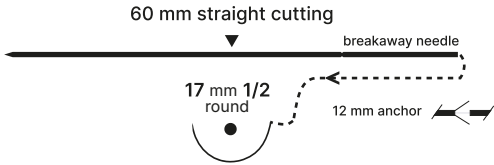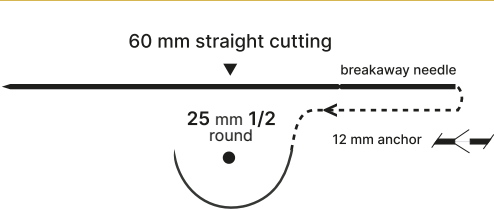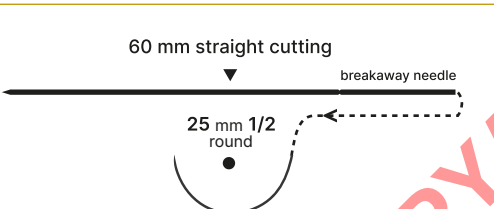60 mm			
straight cutting		60 cm	E5601W
		60 cm	E5603W
		60 cm	E5613W
90 mm			
straight cutting		60 cm	E5903W



Microsurgery & Ophthalmology

Microsurgery

needle information	needle diameter	suture length	0.01	0.1	0.2	0.3	0.4
			12/0	11/0	10/0	9/0	8/0

4 mm								
3/8 round bodied (chord = 3.14 mm)		daylon black	50 μ	13 cm	M0102			
		daylon black	70 μ	13 cm		M1102		
		daylon black	100 μ	13 cm			M2102	M3102
		daylon black	140 μ	13 cm				M4112

5 mm								
3/8 round bodied (chord = 3.92 mm)		daylon black	100 μ	13 cm			M2105	M3105
		propilen blue	100 μ	13 cm			M2205	M3205
		daylon black	140 μ	13 cm				M4115

6 mm								
3/8 round bodied (chord = 4.71 mm)		daylon black	140 μ	13 cm			M2117	M3117
3/8 round bodied double (chord = 4.71 mm)		propilen blue	160 μ	60 cm				P9612

8 mm								
3/8 round bodied, double (chord = 6.27 mm)		propilen blue	160 μ	60 cm				P9812
	 black	propilen blue	160 μ	60 cm				P9814

needle information	needle diameter	suture length	0.2	0.3	0.4
			10/0	9/0	8/0

6 mm							
3/8 spatulated, double (chord = 4.71 mm)		daylon (black)	150 μ	30 cm	M2177	M3167	
		v.silk (blue)	140 μ	30 cm			M4367
		daylon (black)	200 μ	45 cm	M2187		
		v.silk (blue)	200 μ	45 cm			M4387

KONTROLSUZ KOPYADIR . . .

doğsan

IKIH_2026_V0

Doğsan Tıbbi Malzeme Sanayi A.Ş.

Sales and Marketing: Yıldız Cad. No: 55/A Beşiktaş - İSTANBUL

Tel: +90 (212) 258 00 54

E-mail: dogsan@dogsan.com.tr

www.dogsan.com.tr



dogsansurgical



dogsansurgical



DogsanSurgical



dogsansurgical